



Remote license activation and release

FlowLogic uses a floating license, meaning it can be transferred from one computer to another. This is achieved by releasing the license when quitting FlowLogic and activating it, by typing in a username, when launching on a different computer. This process requires the software to connect directly with the license server. However, this may not be possible if network security settings prohibit this communication or if the computer has no Internet connection. To solve this problem, the license can be activated and released remotely.

Please note, remote activation requires a username, which is generated when a license code is registered. If this is the first time you have run FlowLogic, registration can be performed on a different computer that can communicate directly with the license server. Alternatively, please contact Miltenyi tech support for additional assistance.

If network settings prohibit direct communication with the license server but the computer does have access to the Internet via a web browser, then the activation and release can be performed in a few steps.

The general steps to activate a license are:

- 1. Generate a fingerprint unique to your computer**
- 2. Upload the fingerprint to the license activation website and create a key**
- 3. Download the key and import it into FlowLogic**

The general steps to release a license are:

- 1. Generate a fingerprint unique to your computer**
- 2. Upload the fingerprint to the license release website**

In the case where there is no Internet connection, the fingerprint needs to be saved/copied and transferred to a different computer that is online to upload it to the activation and release websites.

Remote Activation – if there is no Internet connection, the Remote Activation window will open automatically. Otherwise, select 'Remote Activation' from the License menu.

1. Type in your username and click generate. A fingerprint is created and displayed.

2. If there is an Internet connection, send the fingerprint directly to the remote activation website.

If there is no Internet connection, copy or save the fingerprint and transfer to a computer that is online and upload the fingerprint to this address.

3. When the fingerprint is submitted, a key will be created. Click to import the key and activate the license.

Remote Activation

FlowLogic could not establish a connection to the license server. To authorize FlowLogic:

- Generate a fingerprint unique to this computer
- Upload the fingerprint to the license server and create a key
- Download the key and import it into FlowLogic

Enter your FlowLogic Username:

Username:

Fingerprint:

1wilsgbk6pum9aw3j8182we8d00kftoqp3nlz6krxbvznzd1-el1uli61zn0s3927fei26togj83rfr6dahya0juuhu4h14xcbf63afz99iw0qllm4749tye727td6rsai4dazy29hi5w9yvwli89unycx3ay49qfmkn26ncbe3wu57jyot0yvgodhz3mzcv17fegspodrhtxjbcwpmgt1s4j2lepjsizv4k944mky5xx61dov3ddbtxqsu6gp0if3tty9oglitnln0fblysvipf9biwe5gr6b

If this computer can access the Internet, after generating a fingerprint, click the following button to automatically download the key via a web browser:

If this computer cannot access the Internet, save the fingerprint to an external drive or network folder. Then, on a computer with Internet access, submit the fingerprint to:
<http://www.inivai.com/check.php?action=activate>


Fingerprint saving options (or copy and paste directly from the window):

Submit the fingerprint to the license key creation website:

If you have generated a FlowLogic key, click the following button:



Remote Activation Website - if uploading automatically from the Remote Activation window, the fingerprint will already be added to the code window. So, simply click 'Submit' to generate a key.

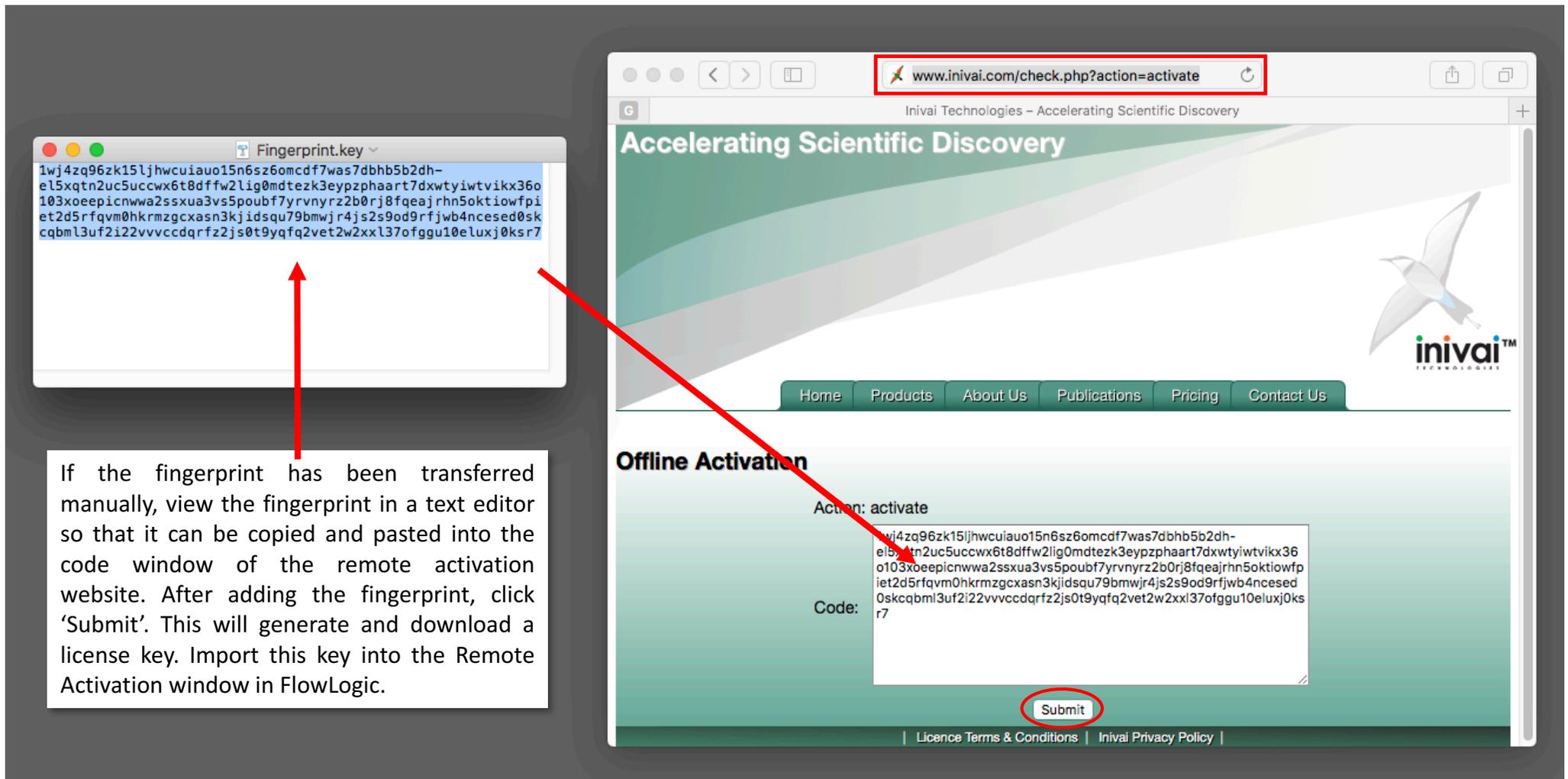

The image shows a screenshot of a web browser displaying the Inivai website. The browser's address bar is highlighted with a red box and contains the URL `www.inivai.com/check.php?action=activate`. The website header features the slogan "Accelerating Scientific Discovery" and the Inivai logo. A navigation menu includes links for Home, Products, About Us, Publications, Pricing, and Contact Us. The main content area is titled "Offline Activation" and contains a form with the following fields:

- Action: activate
- Code: [A text area containing a long alphanumeric fingerprint key]

The "Submit" button at the bottom of the form is circled in red. To the left of the browser window, a separate window titled "Fingerprint.key" displays the same alphanumeric key. A red arrow points from the key in the browser's code field to the key in the text editor window. Another red arrow points from the text editor window to the "Submit" button.


If the fingerprint has been transferred manually, view the fingerprint in a text editor so that it can be copied and pasted into the code window of the remote activation website. After adding the fingerprint, click 'Submit'. This will generate and download a license key. Import this key into the Remote Activation window in FlowLogic.



Remote Release – select 'Remote Release' from the License menu.

A fingerprint is automatically generated. If there is an Internet connection, send the fingerprint directly to the remote release website.

If there is no Internet connection, copy or save the fingerprint and transfer to a computer that is online and upload the fingerprint to this address.

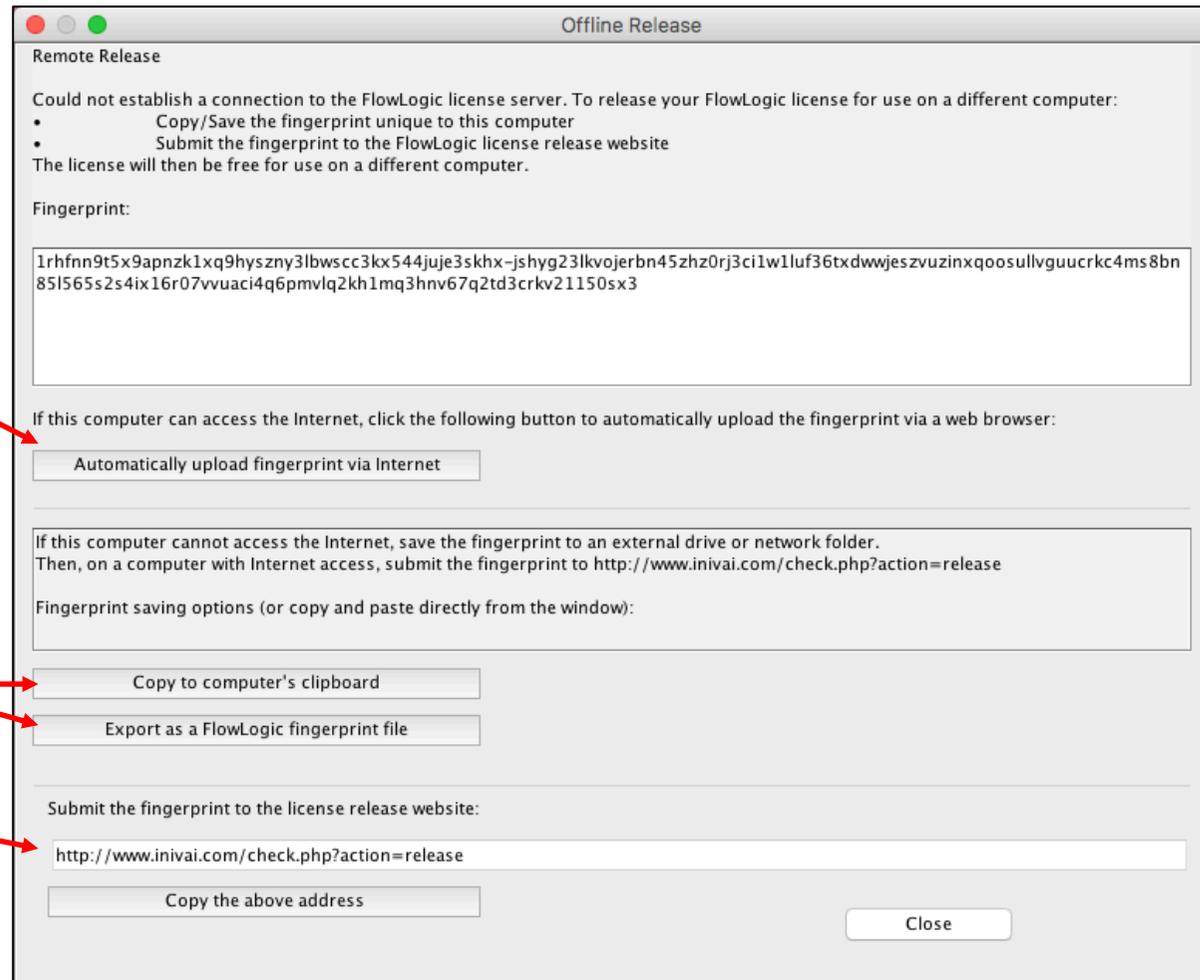

Remote Release Website - if uploading automatically from the Remote Release window, the fingerprint will already be added to the code window. So, simply click 'Submit' to complete the release.

```
Release Fingerprint.key  
1rgptufn3cvzzh3suv3k6t5bv8730tnc8si4npa6ia5pxh39-  
jsauo48kbxytd4oq54q6tp2siprrxpn5p5zvn38kqx323zm0o7xuhio3ik2  
fm49bvjudynotr9g16fvea1ysxx82lmgh76yp98hhepbkhqehguovi22b21  
oia40pww0uwar
```

If the fingerprint has been transferred manually, view the fingerprint in a text editor so that it can be copied and pasted into the code window of the remote release website. After adding the fingerprint, click 'Submit'. This will release the license, making it available for activation in a different computer.

inivai.com/check.php?action=release

Accelerating Scientific Discovery

Home Products About Us Publications Pricing Contact Us

Offline Activation

Action: release

```
1rgptufn3cvzzh3suv3k6t5bv8730tnc8si4npa6ia5pxh39-  
jsauo48kbxytd4oq54q6tp2siprrxpn5p5zvn38kqx323zm0o7xuhio3ik2  
fm49bvjudynotr9g16fvea1ysxx82lmgh76yp98hhepbkhqehguovi22b21  
oia40pww0uwar
```

Code:

Submit

Licence Terms & Conditions | Inivai Privacy Policy

